

## Digital tape recorder



Allows to record sentences (unlimited duration) on 5 separate tracks, for making comparisons of voice quality at different steps of a patient therapy.

---

## Preparation

---

## Equipment

Patient is standing up in front of the microphone.

---

## Software

Launch the SESANE software by clicking this icon in Windows task bar.



In SESANE, enter the patient information :




Patient Informations



Double click on this icon :

## Using Digital Tape Recorder

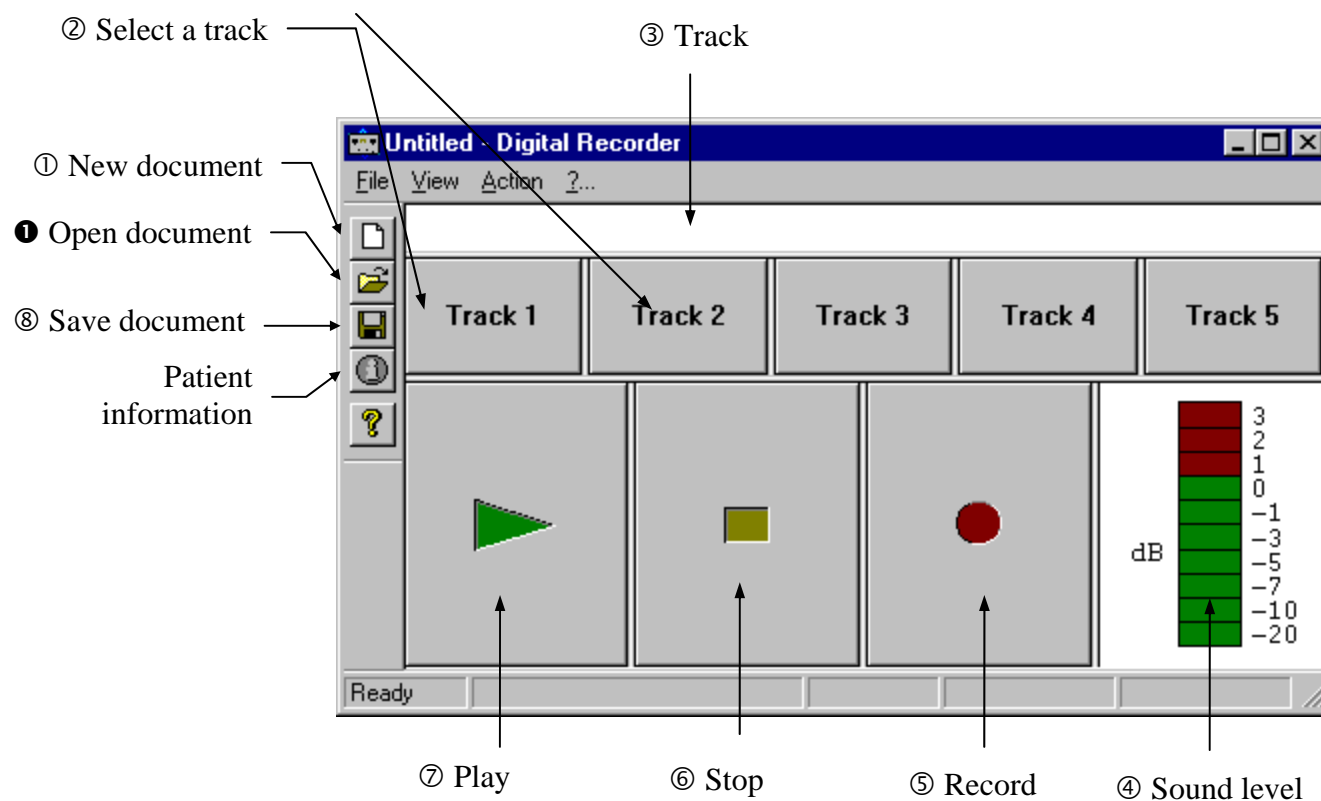

### Procedure

① For a new patient, create a new document

❶ To add a track of a registered patient, open a document.

② Select a track.

③ Take care to select a free track.

④ Speak to test the level. Verify that the sound level is correct (level < 3)

⑤ Click on the Record button. Speak.

⑥ When the patient has finished to speak, click on the Stop button


⑦ Click on the Play button if necessary.

⑧ Save the document

To compare with previous recording :

- select the track ②
- check the track identifier ③
- playback ⑦

## Data Management

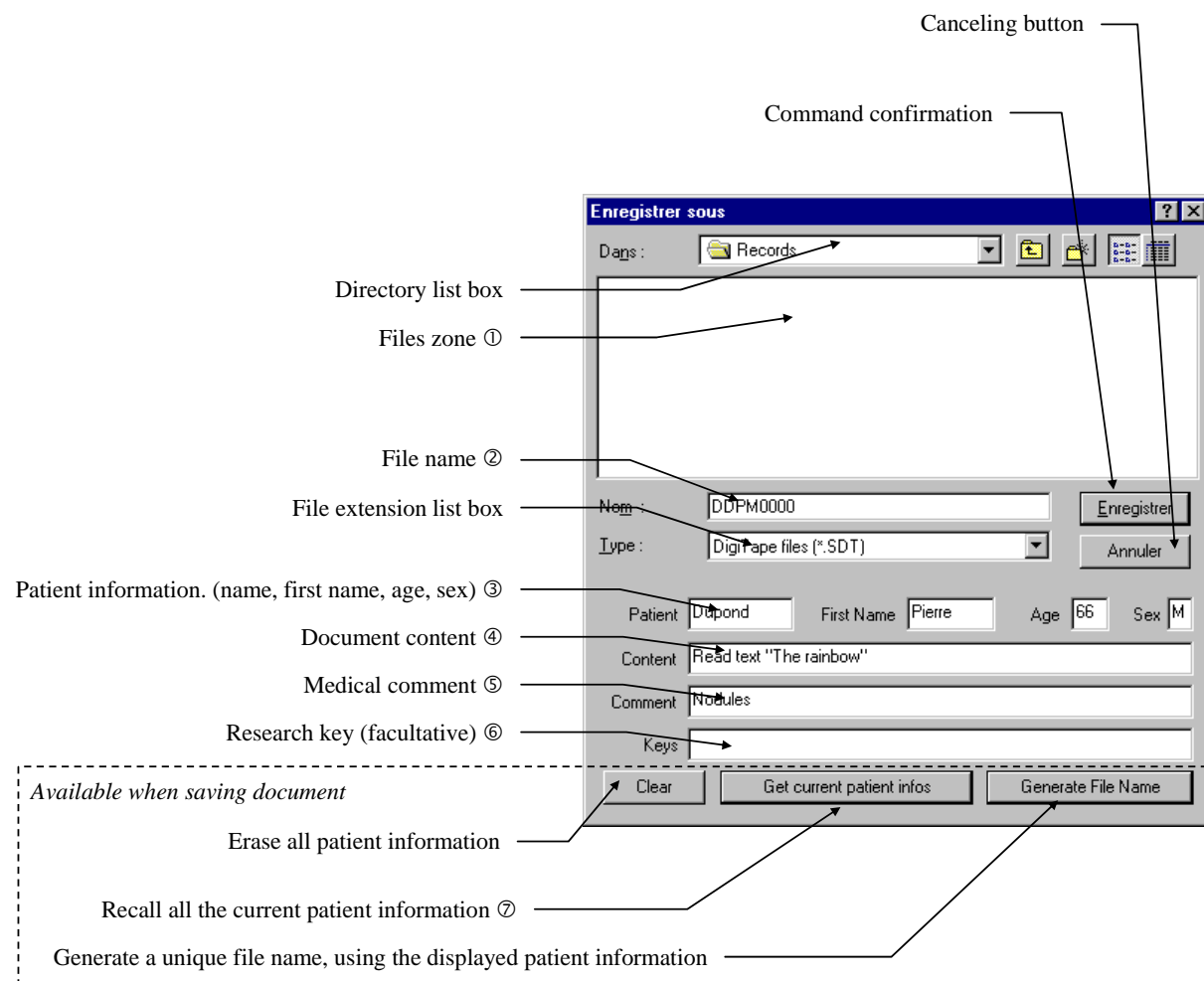

### Save a document

Method 1 : Click on ⑦. The current patient information appear in the fields ③, ④, ⑤, ⑥. A unique filename is automatically generated in ②. Confirm by clicking on ⑨.

Method 2 : Enter manually the patient information in the fields ③, ④, ⑤, ⑥. Click on ⑧. A unique filename appears in ②. Confirm the saving by clicking on ⑨.

Method 3 : Enter manually the patient information in the fields ③, ④, ⑤, ⑥. Enter a file name in ①. Confirm the saving by clicking on ⑨.



### Open a document

Select a document in ① by a single click with the left mouse button.. The file name appears in ② with its information as well in ③, ④, ⑤ ⑥. Confirm your choice by clicking on ⑨.



### To obtain information about the current document

---

---

## Measurement

None.

---

---

## Options

None.